



Miller Keffer Bullock Pedigo

Attorneys and Counselors

110 West Seventh Street
Suite 110
Tulsa OK 74119-1031
918-584-2001
Fax: 918-779-4383

November 7, 2007

Robert George
Kutak Rock LLP
214 W Dickson St
The Three Sisters Building
Fayetteville AR 72701-5221

(VIA EMAIL ONLY)

Re: *State of Oklahoma v. Tyson et al.*,
No. 05-CV-0329-GKF-SAJ

Dear Robert:

You will find attached to this letter Bates stamped copies of the report concerning Identification of a Poultry Specific Biomarker and the November 3, 2007 report entitled Poultry-Specific Biomarker Quantitative PCR Analytical Summary. These are numbered versions of the same reports which I provided to you on Monday November 5, 2007. If you would please substitute these for the earlier production, it might help in keeping these organized.

These documents have been produced in order to assure full compliance with the Court's January 5, 2007 Order (Dkt. 1016). In providing these documents to you, we are not waiving our attorney work product privilege regarding these matters and understand the Court's Order to only relate to the production of facts and not to extend to the production of the opinions or analysis of our experts or to require further invasion of our work product privilege.

As always, should you have any questions concerning this matter, please let me know.

Sincerely,

A handwritten signature in black ink, appearing to read 'L. W. Bullock', is written over a horizontal line.

Louis W. Bullock

Identification of a Poultry Litter Specific Biomarker and Development of a Quantitative Assay

Prepared for:

Dr. Valerie J. Harwood
Dept. Biology, SCA 110
University of South Florida
4202 E. Fowler Ave
Tampa, FL 33620

Prepared by:



North Wind, Inc.
1425 Higham St.
Idaho Falls, Idaho 83402

October 2007

NWI-3517-001

Identification of a Poultry Litter Specific Biomarker and Development of a Quantitative Assay

October 2007

Prepared for:

Dr. Valerie J. Harwood
Dept. Biology, SCA 110
University of South Florida
4202 E. Fowler Ave
Tampa, FL 33620

Prepared by:

North Wind, Inc.
1425 Higham Street
Idaho Falls, Idaho 83402

STOK. 029759

CONTENTS

ACRONYMS	vi
1 INTRODUCTION	1
2 IDENTIFICATION OF POTENTIAL POULTRY LITTER SPECIFIC BIOMARKER	1
2.1 Litter and Soil Community Profiling.....	1
2.1.1 Objective.....	1
2.1.2 Methods	1
2.1.3 Results.....	2
2.2 Cloning and Sequencing of Organisms of Interest.....	3
2.2.1 Objective.....	3
2.2.2 Methods	3
2.2.3 Results.....	3
2.3 Biomarker Sequence Analysis and PCR Primer Design	3
2.3.1 Objective.....	3
2.3.2 Methods	4
2.3.3 Results.....	4
3 TESTING OF THE SPECIFICITY OF THE POTENTIAL POULTRY BIOMARKER TARGETS.....	5
3.1 Test PCR Primers Against Original Soil and Litter Samples	5
3.1.1 Objective.....	5
3.1.2 Methods	5
3.1.3 Results.....	5
3.2 Test PCR Primers Against Closely Related Organisms in BLAST Database.....	7
3.2.1 Objective.....	7
3.2.2 Methods	7
3.2.3 Results.....	7
3.3 Test PCR Primers Against Other Fecal Samples from Within and Outside the Watershed	7
3.3.1 Objective.....	7
3.3.2 Methods	7
3.3.3 Results.....	8
3.4 Sequencing of the Duck and Goose <i>Brevibacteria</i>	11
3.4.1 Objective.....	11
3.4.2 Method	11
3.4.3 Results.....	11

3.5	Test for the Poultry Litter Specific Biomarker in Environmental Samples from Within the Impacted Watershed.....	11
3.5.1	Objective.....	11
3.5.2	Methods	11
3.5.3	Results.....	12
4	OPTIMIZATION AND VALIDATION OF A QUANTITATIVE ASSAY FOR THE POULTRY LITTER SPECIFIC <i>BREVIBACTERIA</i> BIOMARKER	14
4.1	Develop a Quantitative PCR Assay for the <i>Brevibacteria</i> Biomarker	14
4.1.1	Objective.....	14
4.1.2	Method.....	14
4.1.3	Results.....	16
4.2	Test qPCR Specificity to Distinguish Among <i>Brevibacteria</i> spp.	16
4.2.1	Objective.....	16
4.2.2	Methods	17
4.2.3	Results.....	17
4.3	Test for the Poultry Litter Specific Biomarker in Environmental Samples from Within the Impacted Watershed by the qPCR Assay.....	17
4.3.1	Objective.....	17
4.3.2	Method.....	17
4.3.3	Results.....	18
4.4	Determine Dilution Limit of Poultry Litter in Environmental Samples for <i>Brevibacterium</i> Biomarker qPCR Assay	21
4.4.1	Objective.....	21
4.4.2	Methods	21
4.4.3	Results.....	22
5	REFERENCES	23

FIGURES

Figure 1. Outline of MDL protocol.....	15
Figure 2. SYBR green qPCR standard curve of the poultry litter specific primers against plasmid DNA containing the <i>Brevibacteria</i> biomarker DNA. Error bars indicate the standard deviation of cycle thresholds of triplicate samples.....	16
Figure 3. SYBR green melting curve profile of our biomarker sequence (blue) and the <i>Brevibacterium casei</i> [ATCC 35513] (red and green).	17
Figure 4. Outline of "dilution limit" protocol.	21
Figure 5. Copies of <i>Brevibacteria</i> poultry litter biomarker per gram of soil versus grams of soils per litter in the extractions.....	22

TABLES

Table 1. Common T-RFs found in replicate soil and litter samples.	2
Table 2. Primer design for potential biomarkers identified during T-RFLP/clone library assessment.....	4
Table 3. Test of the biomarker-specific PCR analysis of original litter and soil replicates.....	6
Table 4. Results of the PCR amplification of the fecal samples with the four potential poultry litter biomarkers.	9
Table 5. Results of the PCR amplification of the environmental samples for the poultry litter biomarker.....	12
Table 6. Summary of <i>Brevibacteria</i> specific PCR results in DNA extracted from fecal and environmental samples.	14
Table 7. Detection limit of the qPCR assay for a poultry litter specific <i>Brevibacteria</i> in soils and water.	16
Table 8. Results of the test for the poultry litter specific biomarker in environmental samples with the qPCR assay.....	19
Table 9. Results of the PCR, qPCR, and nested qPCR for <i>Brevibacteria</i> in litter.....	22

ACRONYMS

ATCC	American Type Culture Collection
DNA	deoxyribonucleic acid
MDL	method detection limit
MRCF	Molecular Research Core Facility
PCR	polymerase chain reaction
qPCR	quantitative polymerase chain reaction
T-RF	terminal restriction fragment
T-RFLP	terminal restriction fragment length polymorphism
UV	ultraviolet

Identification of a Poultry Litter Specific Biomarker and Development of a Quantitative Assay

1 INTRODUCTION

This report documents the identification of a potential poultry litter specific biomarker, confirmation of the specificity of that biomarker for poultry litter (as compared to other fecal material), and the development of a quantitative polymerase chain reaction (qPCR) assay specific to detection of the poultry biomarker in various environmental media. Methods and results of the initial tests to identify potential poultry litter specific biomarkers are discussed in Section 2. Testing of the sensitivity of the potential biomarkers toward fecally contaminated litter and soil samples, and the specificity of the potential biomarkers against other fecal samples and various environmental media is presented in Section 3. The development and validation of an assay to quantify the poultry litter specific biomarker in various environmental media is discussed in Section 4.

2 IDENTIFICATION OF POTENTIAL POULTRY LITTER SPECIFIC BIOMARKER

The methods utilized to identify a potential poultry litter specific biomarker are discussed in this section. Specifically, the methods and results of microbial community profiling of poultry litter and soils to which the litter had been applied by terminal restriction fragment length polymorphism (T-RFLP) are discussed in Section 2.1. Cloning bacterial DNA from the litter and soil samples and sequencing of the plasmid deoxyribonucleic acid (DNA) containing the terminal restriction fragments (T-RFs) of interest from these clone libraries are presented in Section 2.2. Comparisons of potential biomarker DNA sequences, as compared to published sequences, and development of polymerase chain reaction (PCR) primers specific to these potential biomarkers are discussed in Section 2.3.

2.1 Litter and Soil Community Profiling

2.1.1 Objective

The purpose of the community profiling through the use of T-RFLP was to generate microbial community profiles of all bacteria, *Escherichia coli*, and *Bacteroides* spp. present in the poultry litter and the soils where poultry litter has been applied. The T-RFLP profiles were then used to identify common microorganisms present in both litter and soils on which poultry litter was spread.

2.1.2 Methods

DNA Extraction. DNA was extracted from two poultry litter samples (five replicates each) and two agricultural soil samples (five replicates each) to which the poultry litter was applied. The poultry litter samples from which DNA was extracted were FAC-01A-1 through 5 and FAC-01B-1 through 5, and the soils samples were LAL3-A-2-1 through 5 and LAL3-B-2-1 through 5. Genomic DNA was extracted with Bio101 DNA extraction kits (QBiogene, Inc).

T-RFLP Analysis. T-RFLP was used to generate community profiles targeting *Bacteria* spp., *E. Coli*, and *Bacteroides* spp. The following steps were used during the T-RFLP analysis.

- The extracted DNA was PCR amplified with phosphoramidite fluorochrome 5-carboxyfluorescein (FAM) labeled universal bacterial primers 8F-907R, with *E.coli* genus specific primers (Tsen, et al. 1998), and *Bacteroides* genus specific primers (Bernhard and Field, 2000). All PCR primers targeted the 16S rRNA gene.

- These PCR products were digested with the restriction enzymes *Acil*, *HhaI*, *MspI*, and *HaeIII* (New England Biolabs).
- Common T-RFs for each PCR primer pair (e.g., universal bacteria, *E.coli* or *Bacteroides spp.*), among all 20 poultry litter and soil replicates, were targeted as potential biomarkers of poultry litter.
- T-RFLP analysis was performed by the Idaho State University, Molecular Research Core Facility (MRCF).

2.1.3 Results

Common T-RFs observed in the soil and litter samples are presented in Table 1. No *Bacteroides spp.* fragments were identified as being present in both the soil and litter samples. T-RFs that were well-represented in all samples were selected for further development.

Table 1. Common T-RFs found in replicate soil and litter samples.

<u><i>E. coli</i> PCR products, digested with <i>MspI</i></u>				
T-RF ^a	Litter FAC-01A	Litter FAC-01B	Soil LAL3-A-2	Soil LAL3-B-2
496.0	<u>1,2,4,5</u> ^{b,c}	1,2,3,5	1,2,4	Present in all five
498.9	Present in all five ^d	Present in all five	1,2,4,5	Present in all five
500.8	Present in all five	Present in all five	Present in all five	Present in all five
<u>Universal bacteria PCR products, digested with <i>MspI</i></u>				
80.1	<u>1,2,3,4</u>	Present in all five	Not present in any sample	1,3,4
130.9	1,3,4	Present in all five	3	Not present in any sample
142.9	Present in all five	1,2,3,4	1,4	1
147.3	Present in all five	Present in all five	Present in all five	1,4
158.9	<u>2,3,4</u>	Present in all five	2,3,4,5	1,4
165	1,3,4	Present in all five	1,3,4,5	1,4
a: T-RFs of potential biomarkers are indicated in bold typeface. b: Number indicates the litter or soil replicate sample that the T-RF was identified in. c: An underlined number indicates that the T-RF represented <1% of community in that replicate. d: Indicates that this T-RF was detected in all five subsamples within this sample.				

2.2 Cloning and Sequencing of Organisms of Interest

2.2.1 Objective

The purpose of the cloning and sequencing was to obtain DNA sequences corresponding to the T-RFs of interest found in both poultry litter and soil upon which poultry litter was applied.

2.2.2 Methods

- Clone libraries were constructed from the original genomic DNA extracted from the soil and litter samples and amplified with either universal bacterial or the *E. coli* genus specific primers. Four universal clone libraries were constructed from the following pooled DNA samples (i.e., 1 µl of genomic DNA extract from each sample was added to the PCR reaction for inclusion into the clones): FAC-01A-1 and FAC-01A-4, FAC-01B-3 and FAC-01B-4, LAL3-A-2-1 and LAL3-A-2-4, and LAL3-B-2-1 alone. Four *E. coli* clone libraries were constructed from the following samples: LAL3-A-2-1, FAC-01B-4, FAC-01A-4, and LAL3-B-2-1. The Topo TA cloning kit (Invitrogen) was used for construction of all clone libraries. DNA samples selected for cloning were from those subsamples in which the T-RFs of interest represented a significant portion of the fluorescence profile.
- The plasmids were excised (QIAprep Spin Miniprep Kit, QIAGEN) from the clones and analyzed by T-RFLP (digested with *MspI* alone) to determine which clones contained the T-RFs of interest (see Table 1).
- Plasmids containing the T-RFs of interest were amplified by PCR and sequenced using the primers T7, T3, 519R, and 338F for double coverage of the 16s rDNA.

2.2.3 Results

- Only three of the six T-RFs representing potential biomarkers were found in the universal clone library (i.e., T-RFs 142.9, 147.3, and 158.9). Sequences representing all three *E. coli* T-RFs were found in the *E. coli* clone library.
- After developing two clone libraries and screening an additional 88 clones with T-RFLP, the target biomarkers with T-RFs of 80.1, 130.6 and 165 were not found. A total of 350 clones were screened.
- T-RFLP sequence analysis and DNA sequencing of each clone was performed by the MRCF.

2.3 Biomarker Sequence Analysis and PCR Primer Design

2.3.1 Objective

The purpose of the biomarker sequence analysis was to compare the potential biomarker DNA sequences to published 16S rRNA DNA sequences to determine whether these sequences had previously been observed, and whether they were identified as a particular species or associated with a particular host animal. Additionally, these sequences were used to design PCR primers specific to each biomarker.

2.3.2 Methods

- The three universal and the three *E. coli* sequences corresponding to the T-RFs of interest were compared to the BLAST database (i.e., National Center for Biotechnology Information) to determine closely related organisms to our potential biomarkers and sites amenable for the design of PCR primers.
- PCR primers were designed for the three biomarkers from the universal bacterial library and one of the *E. coli* biomarkers, targeting regions of variability between our sequence and the database sequences of the top 20 matches in the BLAST database. PCR primers were designed using the Primer Express v2.0 software (Applied Biosystems). These primers were analyzed for thermodynamic folding problems and compared to the RDP-II database (Michigan State University) to determine what other organisms they might amplify. The results of the analysis of the forward and reverse PCR primer comparison against the RDP-II database are presented in the Table 2.

Table 2. Primer design for potential biomarkers identified during T-RFLP/clone library assessment.

Potential Biomarker Clone Number/ T-RF/Organism *	Forward Primer Accession Number of Closest Match in Sequence Database RPDII	Reverse Primer Accession Number of Closest Match in Sequence Database RPDII
Clone SA19 T-RF 158.9 <i>Kineococcus</i>	AY919955, AY426452, AF195447, AF513961, AY862810, AM085954, AM182287, AM182297, AM182298, AM182299	Primer sequence did not match any organisms in the database
Clone SB37 T-RF 142.9 Uncultured organisms	Primer sequence did not match any organisms in the database	Primer sequence did not match any organisms in the database
Clone LA35 T-RF 147.3 <i>Brevibacterium</i>	Primer sequence did not match any organisms in the database	Primer sequence did not match any organisms in the database
Clone SA15 T-RF 500.8 <i>Pantoea ananatis</i>	Primer sequence did not match any organisms in the database	AJ010486, DQ221344, AF364845, AF364844, AY528223, AY579209, AY579211, U80196, U80209, AB004758, AB027693, AY530796, AJ629190, AB178169, AB178170, AY898643, AB114622, DQ133548, DQ195524, AB242937, AB242945, AB242946, AB242979 And an additional 30 sequences <i>Enterobacter</i> spp.
* organisms were identified based on a BLAST search		

2.3.3 Results

- Only one *E. coli* biomarker had a variable region observed in all 20 matches of the closest related organisms in the BLAST database. PCR primers were designed for this variable

region and another region that had a mismatch between our sequence and the database sequence in 10 of the top 20 matching sequences.

- One primer pair was designed that was determined to be specific for T-RF of interest 147.3 from clone LA-35, a *Brevibacterium* spp., and would have no matches compared to the sequences in the RDPII database.
- One primer pair was designed for a *Kineococcus* spp. corresponding to T-RF 158.9, where the forward primer matched the sequences of 10 other organisms in the RDPII database, but the reverse primer was specific to this *Kineococcus* spp.
- One primer pair was designed for an organism not matching any cultured organisms in the BLAST or RDPII databases corresponding to T-RF 142.9.

3 TESTING OF THE SENSITIVITY AND SPECIFICITY OF THE POTENTIAL POULTRY BIOMARKER TARGETS

The methods utilized to test the sensitivity and specificity of the potential poultry litter specific biomarkers are discussed in this section. Specifically, the methods and results for the testing the PCR primers in the original soil and litter samples used for the T-RFLP are presented in Section 3.1, while the testing of the PCR primers for the biomarkers against closely related organisms in the BLAST database search are presented in Section 3.2. Testing of the PCR primers for the various potential biomarkers against other fecal material is presented in Section 3.3. Cloning and sequencing of PCR amplicons derived from fecal samples that amplified with the *Brevibacteria* specific primers and comparison to the biomarker sequence is presented in Section 3.4 Finally, testing of the PCR primers of the confirmed poultry litter specific biomarker in environmental samples from the potentially poultry litter impacted watershed are presented in Section 3.5.

3.1 Test PCR Primers Against Original Soil and Litter Samples

3.1.1 Objective

The objective of this test was to determine if the PCR primers specific for the various biomarkers amplified DNA from the original contaminated soil and poultry litter samples used to find the T-RFs of interest (i.e., potential biomarker sequences).

3.1.2 Methods

- A nested PCR approach was used to increase sensitivity of the PCR assay by first amplifying the DNA with the universal bacterial primers 8F-907R or *E. coli* species-specific primers, purifying the PCR products (QIAquick PCR purification kits, QIAGEN), and then amplifying the DNA with the biomarker-specific primers.
- PCR products were separated by gel electrophoresis, stained with ethidium bromide, and visualized by ultraviolet (UV) light.

3.1.3 Results

The results of the nested PCR with the potential biomarker-specific PCR primers of DNA from the original litter and soil samples are presented in Table 3.

Table 3. Test of the biomarker-specific PCR analysis on original litter and soil replicates.

	Clone LA35 <i>Brevibacterium</i> spp.		Clone SB37 Unknown genus		Clone SA19 <i>Kineococcus</i> spp.		Clone SA15 <i>E. coli</i>	
Sample	% of T-RF Profile	Amplified with LA35 Primers?	% of T-RF Profile	Amplified with SB37 Primers?	% of T-RF Profile	Amplified with SA19 Primers?	% of T-RF Profile	Amplified with SA15 Primers?
Litter Sample Results								
FAC-01-A-1	2.6	yes	2.8	yes	not present ^b	not run	26.4	yes
FAC-01-A-2	2.6	yes	4.6	no	0.8	yes	13	no
FAC-01-A-3	3.4	yes	5.4	yes	1	not run ^c	no data ^a	yes
FAC-01-A-4	3.3	yes	5.9	yes	1.3	not run ^c	23.3	no
FAC-01-A-5	no data ^a	yes	no data ^a	not run ^c	no data ^a	not run ^c	28.5	yes
FAC-01-B-1	3.3	yes	3.6	not run ^c	1.5	not run ^c	33	yes
FAC-01-B-2	4	yes	5.5	yes	1.1	no	43.8	yes
FAC-01-B-3	3.5	yes	6.8	yes	1.3	yes	16.5	yes
FAC-01-B-4	3.2	yes	5.7	not run ^c	1.2	not run ^c	29.1	no
FAC-01-B-5	4.5	yes	8	not run ^c	1.1	not run ^c	48.5	no
Soil Sample Results								
LAL3-A-2-1	7.1	yes	0.8	yes	2.3	not run ^c	6.7	yes
LAL3-A-2-2	12.7	no	not present ^b	yes	3.4	not run ^c	14.3	no
LAL3-A-2-3	9	yes	1	yes	3.6	not run ^c	25.3	yes
LAL3-A-2-4	6.9	yes	0.8	yes	3.2	yes	10.4	yes
LAL3-A-2-5	9.5	no	0.9	yes	3.6	not run ^c	2.5	yes
LAL3-B-2-1	6	yes	0.8	yes	3	yes	18.1	yes
LAL3-B-2-2	no data ^a	yes	no data ^a	yes	no data ^a	yes	6.8	yes
LAL3-B-2-3	not present ^b	yes	not present ^b	yes	3	no	9.1	yes
LAL3-B-2-4	6.3	yes	0.8	yes	3.5	no	2.6	yes
LAL3-B-2-5	no data ^a	yes	no data ^a	yes	no data ^a	yes	7.1	yes

	Clone LA35 <i>Brevibacterium</i> spp.		Clone SB37 Unknown genus		Clone SA19 <i>Kineococcus</i> spp.		Clone SA15 <i>E. coli</i>	
Sample	% of T-RF Profile	Amplified with LA35 Primers?	% of T-RF Profile	Amplified with SB37 Primers?	% of T-RF Profile	Amplified with SA19 Primers?	% of T-RF Profile	Amplified with SA15 Primers?
<p>a <i>No data</i> indicates that the T-RFLP analysis was not completed on this sample.</p> <p>b <i>Not present</i> indicates that the organism represented by that T-RF was not found in the original analysis (see Section 2.1).</p> <p>c <i>Not run</i> indicates that this sample was not run with PCR.</p>								

3.2 Test PCR Primer Set LA35 Against a Closely Related Bacterium

3.2.1 Objective

The objective of this test was to determine if the PCR primers for the *Brevibacterium* LA35 potential poultry litter biomarker amplified the same product in *Brevibacterium* sp. CHNDP32 (DQ337537), the fourth closest related organism found in the BLAST search.

3.2.2 Methods

- 16S rDNA of the closely related organism identified in the BLAST search was obtained from Dr. Chee-Sanford from the University of Illinois at Urbana-Champaign.
- The *Brevibacterium* sp. CHNDP32 was PCR amplified using the LA35 primers, and PCR products were separated by gel electrophoresis, stained with ethidium bromide, and visualized by UV light.

3.2.3 Results

The *Brevibacterium* biomarkers primers did not amplify the *Brevibacterium* sp. CHNDP32.

3.3 Test PCR Primers Against Other Fecal Samples from Within and Outside the Watershed

3.3.1 Objective

The purpose of this test was to determine the specificity of the four potential poultry biomarker targets (LA35, SB37, SA19 and SA15) against other sources of fecal contamination from within and outside the potentially poultry impacted watershed.

3.3.2 Methods

- Fecal samples were collected in duplicate from beef and dairy cattle, swine, geese, ducks, and humans from inside and outside the potentially poultry impacted watershed. Field blank controls were included with each type of fecal sample. The fecal samples were preserved in glycerol and shipped on ice to the laboratory.
- Genomic DNA was extracted from all fecal samples using the Bio 101 FastDNA SPIN Kit for Soil, DNA extraction kit (Qbiogene, Inc).

- Samples were tested for PCR amplification with the four potential biomarkers, and products were separated by gel electrophoresis, stained with ethidium bromide, and visualized by UV light. Samples were run with a nested PCR approach.

3.3.3 Results

The results of the PCR amplification of the fecal samples with the four potential poultry litter biomarkers are presented in Table 4.

Table 4. Results of the PCR amplification of the fecal samples with the four potential poultry litter biomarkers.

Sample	Type of Fecal Sample	Inside or Outside the Watershed	City	Did the Potential Biomarker Sequence PCR Amplify in the Fecal Sample?			
				Clone LA35 <i>Brevibacterium</i>	Clone SB37 Unknown Organism	Clone SA19 <i>Kleococcus</i>	Clone SA15 <i>E. coli</i>
MAN-BC-1-a	Beef Cattle	Outside		No	No	No	No
MAN-BC-1-b		Outside		No	No	No	No
MAN-BC-2-a		Outside		No	Yes	No	No
MAN-BC-2-b		Outside		No	Yes	No	No
MAN-BC-3-a		Outside		No	No	No	No
MAN-BC-3-b		Outside		No	No	No	No
MAN-BC-4-a		Outside		No	No	No	No
MAN-BC-4-b		Outside		No	No	No	No
MAN-BC-5-a		Outside		No	No	No	No
MAN-BC-5-b		Outside		No	Yes	Yes	No
MAN-BC-6-a		Outside		No	No	No	No
MAN-BC-6-b		Inside		No	Yes, faint ^b	No	No
MAN-BC-7-a		Inside		No	No	No	No
MAN-BC-7-b		Inside		No	Yes	Yes, faint ^b	Yes
MAN-BC-8-a		Inside		No ^a	Yes	Yes	Yes
MAN-BC-8-b		Inside		No	No	Yes	No
MAN-BC-9-a		Inside		No	No	Yes, faint ^b	No
MAN-BC-9-b		Inside		No	No	No	No
MAN-BC-10-a		Inside		No	Yes, faint ^b	Yes, faint ^b	No
MAN-BC-10-b		Inside		No	Yes	Yes	No
MAN-BC-F-a	Field Blank	Outside		No	Yes	No	No
MAN-DC-1	Dairy Cattle	Outside		No	Yes	Yes, faint ^b	Yes
MAN-DC-2-a		Outside		No	No	No	No
MAN-DC-2-b		Outside		No	No	No	No
MAN-DC-3		Inside		No	Yes	No	No
MAN-DC-3-b		Inside		No	Yes	No	No
MAN-DC-F	Field Blank	Outside		No	No	No	No
MAN-SW-1-a	Swine	Outside		No	Yes	Yes	No
MAN-SW-1-b		Outside		No	Yes	No	No
MAN-SW-2		Inside		No	No	No	No

Poultry Litter Biomarker
NWL-3517-001

9

North Wind, Inc.
October 2007

Table 4. (continued).

Sample	Type of Fecal Sample	Inside or Outside the Watershed	City	Did the Potential Biomarker Sequence PCR Amplify in the Fecal Sample?			
				Clone LA35 <i>Brevibacterium</i>	Clone SB37 Unknown Organism	Clone SA19 <i>Kleococcus</i>	Clone SA15 <i>E. coli</i>
MAN-DK-1-a	Duck	Outside		No	Yes	Yes	Yes
MAN-DK-1-b		Outside		Yes ^a	Yes	Yes	Yes
MAN-DK-2-a		Outside		No	Yes	Yes	Yes
MAN-DK-2-b		Outside		No	Yes	Yes	Yes
MAN-DK-3-a		Inside		No	Yes	Yes	Yes
MAN-DK-3-b		Inside		No	No	Yes	Yes
MAN-DK-4-a		Inside		No	No	No	No
MAN-DK-4-b		Inside		No	No	No	No
MAN-DK-5-a		Inside		No	Yes	No	No
MAN-DK-5-b		Inside		No	No	No	Yes
MAN-DK-F	Field Blank	Inside		No	No	No	No
MAN-GS-1-a	Goose	Outside		Yes, faint ^{a,b}	Yes	Yes	Yes
MAN-GS-1-b		Outside		No	Yes	Yes, faint ^b	No
MAN-GS-2-a		Outside		No	Yes	No	No
MAN-GS-2-b		Outside		No	No	No	No
MAN-GS-3-a		Inside		No	No	No	No
MAN-GS-3-b		Inside		No	Yes	Yes	Yes
MAN-GS-4-a		Inside		No	Yes	No	No
MAN-GS-4-b		Inside		No	Yes	No	No
MAN-GS-5-a		Outside		No	Yes	Yes	Yes
MAN-GS-5-b		Outside		No	Yes	Yes	Yes
MAN-HM-1	Waste Water Treatment Plant	Outside	Claremore	No	Yes	Yes	Yes
MAN-HM-2		Inside	Siloam Springs	No	Yes	Yes	Yes
MAN-HM-3		Inside	Fayetteville	No	Yes	No	No
MAN-HM-4	Septic System	Outside	Tulsa	No	Yes	Yes	No
MAN-HM-5		Inside	Fayetteville	No	Yes	No	No
MAN-HM-6		Inside	Siloam Springs	No	No	No	No

a: Samples were re-extracted and re-run in duplicate to confirm result.
b: Faint, indicates a very weak band was observed on the gel.

Poultry Litter Biomarker
NWL-3517-001

3.4 Sequencing of the Duck and Goose Amplicon Derived from the LA35 Primer Set

3.4.1 Objective

The purpose of the cloning and sequencing of the duck and goose amplicon derived from the LA35 primer set was to determine if these DNA sequences contained variable regions that could be targeted for a more specific set of PCR primers.

3.4.2 Method

- Several clone libraries were constructed from the duck sample MAN-DK-1-b and the goose sample MAN-GS-1-a by PCR amplifying DNA from these fecal samples with the LA35 primer set, with primers 8F-907R and with 8F-1492R, and cloning with the TOPO-TA cloning kit (Invitrogen). The plasmids were excised from these clones (QIAprep Spin Miniprep Kit, QIAGEN).
- Extracted plasmids were DNA sequenced using the PCR primers T7, T3, 8F, 907R, or 1492R (as appropriate) by the MRCF. Sequences were compared to the poultry litter *Brevibacterium* biomarker using BioEdit V.7.0.5.3 to look for variable regions between the DNA sequences.

3.4.3 Results

The regions of the 16S rRNA gene of the duck and goose DNA targeted by the poultry litter specific *Brevibacterium* biomarker PCR primers are identical to our biomarker *Brevibacterium* sequence.

After screening 20 clones from each of the two clone libraries developed (a total of 40 clones) from the duck and the goose fecal samples, only one organism was identified as a *Brevibacterium* using the BLAST database. This sequence was identical to the *Brevibacterium* biomarker sequence obtained from the original soil and litter samples within the 1,200 DNA base pairs sequenced.

3.5 Test for the Poultry Litter Specific Biomarker in Environmental Samples from Within the Impacted Watershed

3.5.1 Objective

The purpose of this test was to determine if the *Brevibacterium* biomarker that is specific to poultry litter could be detected in environmental media (poultry litter, soil and water samples) from within a potentially poultry litter impacted watershed.

3.5.2 Methods

- Environmental samples were collected from within the potentially impacted watershed and included chicken and turkey litter samples, soil samples from which the litter had been applied, edge of field runoff water samples from the fields to which the litter had been applied, river water to which the runoff samples drained, and lakes downgradient from the previously collected river samples. Additionally, groundwater samples within the potentially impacted watershed were collected and analyzed.

- Water samples were collected in duplicate 1-L sterile nalgene bottles and shipped on ice to the laboratory, where they were immediately filtered. The filters containing the microorganisms were frozen at -80°C until DNA extraction.
- Genomic DNA was extracted from the filters and soil and litter samples using the Bio 101 FastDNA SPIN Kit for Soil, DNA extraction kit (Qbiogene, Inc).
- Samples were tested for PCR amplification with the LA35 PCR primers, and products were separated by gel electrophoresis, stained with ethidium bromide, and visualized by UV light. Samples were run with a nested PCR approach

3.5.3 Results

The results of the PCR amplification of the environmental samples for the poultry litter biomarker are presented in Table 5.

Table 5. Results of the PCR amplification of the environmental samples for the poultry litter biomarker.

Sample	Type of Sample	Amplified with LA35 PCR Primers?
FAC1-6-20-06	Litter	Yes
FAC2-6-21-06	Litter	Yes
FAC-03-7-6-06	Litter	Yes
FAC-4-7-12-06	Litter	Yes
FAC-5-7-13-06	Litter	Yes
FAC-6-7-20-06	Litter	No ^a
FAC-7-8-3-06	Litter	Yes
FAC-8-8-15-06	Litter	Yes
FAC-9-8-31-06	Litter	Yes
LAL5-A-2-6-13-06	Soil	No
LAL5-C-2-6-12-06	Soil	Yes
LAL-7-A-2-6-20-06	Soil	Yes
LAL-7-B-2-6-20-06	Soil	Yes, faint ^b
LAL-7-C-2-6-19-06	Soil	Yes, faint ^b
LAL10-B-2-6-26-06	Soil	No
LAL10-A-2-6-26-06	Soil	No
LAL10-A-4-6-26-06	Soil	No
LAL8-A-2-6-19-06	Soil	Yes
LAL8-B-2-6-21-06	Soil	No
LAL9-D-2-6-22-06	Soil	No
LAL9-B-2-6-22-06	Soil	Yes
LAL9-A-2-6-22-06	Soil	Yes
LAL7-D-2-6-29-06	Soil	Yes

Table 5. (continued).

Sample	Type of Sample	Amplified with Poultry Litter Specific <i>Brevibacteria</i> PCR Primers?
LAL8-D-2-6-20-06	Soil	Yes
LAL11-C-2-6-28-06	Soil	No
LAL11-D-2-6-28-06	Soil	No
LAL11-A-2-6-29-06	Soil	Yes
LAL11-D-2-Q-6-28-06	Soil	Yes
LAL12-A-2-7-6-06	Soil	No
LAL12-A-2-Q-7-6-06	Soil	Yes
LAL12-C-2-7-7-06	Soil	No
LAL12-D-2-7-7-06	Soil	No
LAL 13-A-2-7-6-06	Soil	Yes
LAL 13-C-2-7-7-06	Soil	No
LAL 13-C-2-Q-7-7-06	Soil	Yes
LAL 13-D-2-7-6-06	Soil	No
EOF-1-6-17-06	Water	No
EOF-Q2-6-17-06	Water	No
EOF-Q1-6-17-06	Water	Yes
EOF-SPREAD073B-6-18-06	Water	Yes
EOF-SPREAD023-6-18-06	Water	Yes
EOF-SPREAD044-6-18-06	Water	No
EOF-SPREAD068-6-18-06	Water	Yes
a: This sample contained a high percentage of soils and very little "litter." b: Faint indicates that a weak band was visible on the gel.		

A summary of the presence or absence of the *Brevibacteria* biomarker in all samples analyzed to date is included in Table 6.

Table 6. Summary of *Brevibacteria* specific PCR results in DNA extracted from fecal and environmental samples.

Type of Sample	Number Analyzed	Number Positive	% Containing LA35 sequence	Note
Original soil and litter samples	20	18	90.0	
Fecal Samples	57	2	3.5	LA35 sequences were present in one duck and one goose sample from outside the watershed
Additional litter samples	9	8	88.9	One of the litter samples (FAC-06) was taken from a poultry house without a concrete floor and likely contained a high percentage of soil.
Additional soil samples	27	14	51.9	
Edge of field water samples	7	4	57.1	

4 OPTIMIZATION AND VALIDATION OF A QUANTITIVE ASSAY FOR THE POULTRY LITTER SPECIFIC *BREVIBACTERIA* BIOMARKER

4.1 Develop a Quantitative PCR Assay for the *Brevibacterium* LA35 Biomarker

4.1.1 Objective

The purpose of this work was to develop a qPCR assay for quantification of the LA35 biomarker in various environmental and fecal samples.

4.1.2 Method

- LA35 PCR primers developed in Section 2 and validated for specificity to poultry litter in Section 3 were applied as qPCR primers using SYBR Green chemistry (dye) on a Chromo4 qPCR system (Bio-Rad).
- Plasmids containing the LA35 DNA sequence were used to create a standard curve, to determine a method detection limit (MDL) for positive control plasmids, and to determine the efficiency of the reaction. All standard curves were run in triplicate to verify the reproducibility of assay.
- Detection limits of the assay in environmental samples were determined by spiking soil and water samples with LA35 contained on plasmids. Controls were run to determine the average number of plasmids present in the *E. coli* cells used as plasmid carriers. Additionally, nanopure water and a composite water sample from the watershed were spiked with the *E. coli* containing the biomarker sequence and were filtered according to

the standard filtering methods. Cloning was performed using the Topo TA cloning kit (Invitrogen). Genomic DNA was extracted from the spiked samples using the Bio 101 FastDNA SPIN Kit for Soil, DNA extraction kit (Qbiogene, Inc). A graphic depicting the MDL protocol is presented in Figure 1.

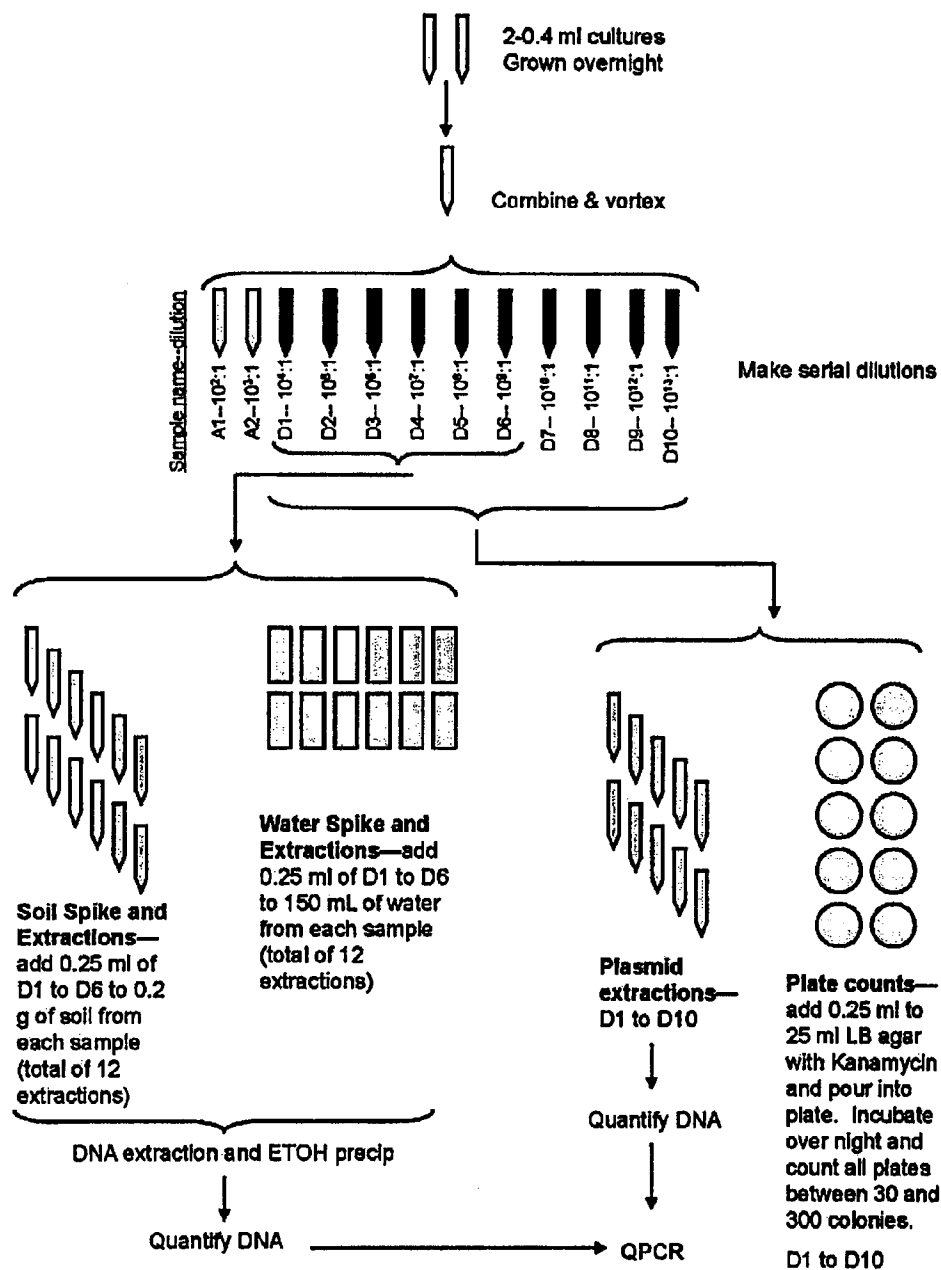


Figure 1. Outline of MDL protocol.

4.1.3 Results

The SYBR green qPCR standard curve is presented in Figure 2. Additionally, the plasmid MDL and reaction efficiency are presented in the figure. Efficiency of the qPCR reaction was determined by equation 1.

$$\text{Efficiency} = -1 + 10^{(1-1/\text{slope})}$$

(1)

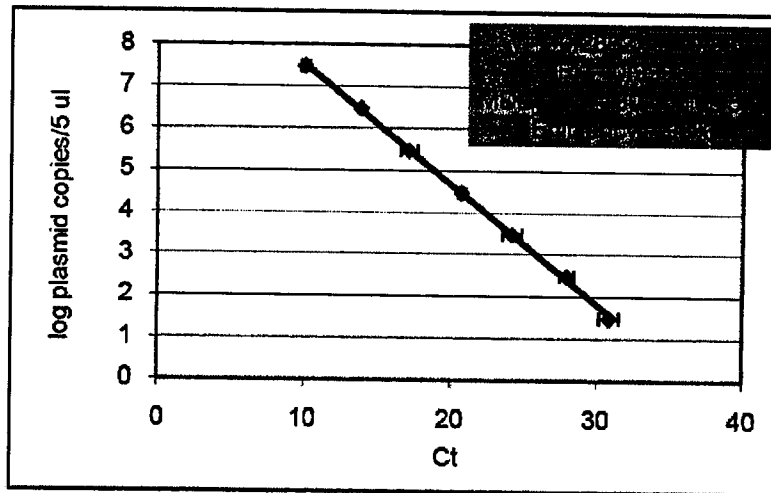


Figure 2. SYBR green qPCR standard curve of the poultry litter specific primers against plasmid DNA containing the *Brevibacteria* biomarker DNA. Error bars indicate the standard deviation of cycle thresholds of triplicate samples.

The results of the minimum detection limit test in spiked environmental samples are presented in Table 7.

Table 7. Detection limit of the qPCR assay for a poultry litter specific *Brevibacteria* in soils and water.

Sample Type	Minimum Detection Limit	Units
Plasmid DNA (standard)	6	copies/ μ L DNA extraction
Nanopure Water	18	cells/L
Composite Water Sample	78	cells/L
LAL11D-2Q-6-28-06 (soil) with sepharose cleanup	73	cells/g

4.2 Test qPCR Specificity to Distinguish Among *Brevibacteria* spp.

4.2.1 Objective

The purpose of this test was to determine if the qPCR assay is specific enough to distinguish between the poultry litter biomarker and the closely related (but not identical) *Brevibacterium casei* 16S rRNA gene.

4.2.2 Methods

A *Brevibacterium casei* culture was ordered from American Type Culture Collection (ATCC) and the DNA was extracted with the standard protocol. The extracted DNA was then tested for amplification with our qPCR protocol.

4.2.3 Results

Our qPCR primers amplified the *Brevibacterium casei* 16S rRNA gene. As shown in Figure 3, we are able to distinguish *B. casei* from the LA35 sequence by the SYBR green melt curves.

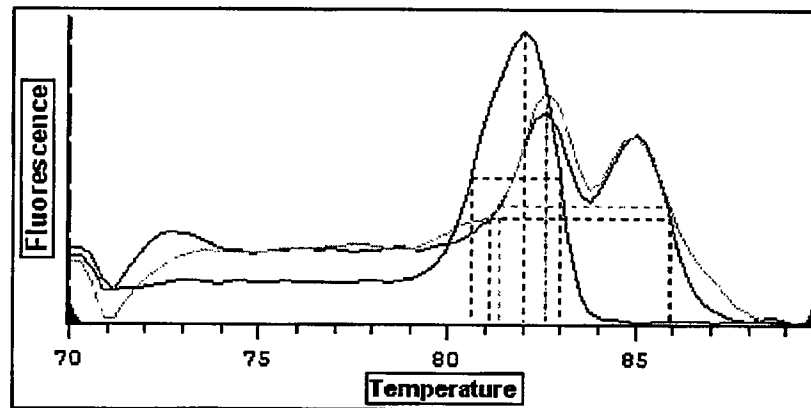


Figure 3. SYBR green melting curve profile of the LA35 16S rRNA sequence (blue) and that of *Brevibacterium casei* [ATCC 35513] (red and green).

4.3 Test for the Poultry Litter Specific Biomarker in Environmental Samples from Within the Impacted Watershed by the qPCR Assay

4.3.1 Objective

The purpose of this test was to determine if the qPCR assay could detect the poultry litter specific biomarker in environmental samples from within the potentially impacted watershed.

4.3.2 Method

- Genomic DNA was extracted from various environmental media (poultry litter, soil and water samples) per the previously described methods. Samples were selected that contained varying levels of fecal indicator bacteria, which was used to gauge the expected biomarker concentration (Table 8)
- DNA was subjected to a diagnostic PCR to verify that the samples did not contain contaminants that might inhibit qPCR and to determine the appropriate sample volume to use for the qPCR assay.
- DNA was then analyzed by the qPCR assay for the poultry litter specific biomarker. Any sample not amplifying in the qPCR assay was tested by a nested qPCR assay, where

universal bacterial primers were first used to amplify the 16S rRNA gene, and then these PCR products were analyzed by the qPCR assay.

4.3.3 Results

The results of the qPCR and nested qPCR assay are presented in Table 8.

Table 8. Results of the test for the LA35 biomarker in environmental samples with the qPCR assay.

Sample ID	Matrix	Expected Biomarker Concentration	DNA (ng/L or ng/g)*	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter)		qPCR Matrix Spike Amplified?*	Nested qPCR Amplified?*	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
				\pm	\pm				
EOF-spr-010-5-9-06	Water	High	1.7	1.05E+07	\pm	1.70E+06	Yes	Yes	No
EOF-spr-17A-01-5-1-06	Water	High	72.5	2.48E+06	\pm	4.71E+05	Yes	Yes	Yes
EOF-spr-023-6-18-06	Water	High	4.3	1.11E+05	\pm	2.49E+03	Yes	Yes	No
EOF-spread-073B-6-18-06	Water	High	19.2	1.92E+06	\pm	4.42E+04	Yes	Yes	No
LAL16-SPR2-7-18-06	Water	High	-1.0	Not present			Yes	N/A	N/A
LAL16C-2-7-18-06	Soil	High	9.5	1.42E+04	\pm	1.97E+03	Yes	Yes	No
LAL11C-2-6-28-06	Soil	High	73.2	Present, not quantifiable			Yes	Yes	No
HFS16-BF1-01-6-15-06	Water	Medium	6.8	4.00E+03	\pm	1.60E+03	Yes	Yes	No
SALspr-6-28-06	Water	Medium	-0.6	5.82E+02	\pm	1.56E+02	Yes	Yes	No
LAL15-SP2-7-11-06	Water	Medium	5.0	2.89E+03	\pm	7.69E+02	Yes	Yes	No
RS-PRICEck-01-4-29-06	Water	Medium	4.7	3.45E+05	\pm	1.43E+05	Yes	Yes	No
RS-574-BIO	Water	Medium	6.7	1.80E+05	\pm	6.09E+04	Yes	Yes	No
Lk04-0-01-5-16-06	Water	Low	6.8	3.69E+03	\pm	3.24E+03	Yes	Yes	No
HFS28A-BF1-01-6-15-06	Water	Low	-0.7	2.48E+03	\pm	1.28E+03	Yes	Yes	Yes
Rs-1-01-8-8-06	Water	Low	7.0	3.19E+04	\pm	6.75E+03	Yes	Yes	Yes
FAC-01A-1	Litter	High	33.7	2.18E+09	\pm	3.53E+08	Yes	Yes	No
FAC-01A-2	Litter	High	4.7	2.47E+08	\pm	3.22E+07	Yes	Yes	No
FAC-01A-3	Litter	High	-0.5	2.67E+07	\pm	2.69E+06	Yes	Yes	No

Poultry Litter Biomarker
NWT-3517-001

19

North Wind, Inc.
October 2007

Table 8. (continued).

Sample ID	Matrix	Expected Biomarker Concentration	DNA (ng/L or ng/g)	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter)		qPCR Matrix Spike Amplified?*	Nested qPCR Amplified?*	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
FAC-01A-4	Litter	High	3.4	1.49E+08	\pm	1.10E+07	Yes	Yes	No
FAC-01A-5	Litter	High	4.1	5.67E+08	\pm	3.75E+07	Yes	Yes	No
FAC-01B-1	Litter	High	94.5	3.94E+09	\pm	6.28E+08	Yes	Yes	No
FAC-01B-2	Litter	High	40.5	2.66E+09	\pm	7.57E+08	Yes	Yes	No
FAC-01B-3	Litter	High	34.5	4.75E+06	\pm	4.23E+06	Yes	Yes	No
FAC-01B-4	Litter	High	117.1	5.99E+09	\pm	1.74E+09	Yes	Yes	No
LAL8-A-2-6-19-06	Soil	High	22.34	7.00E+03	\pm	4.43E+02	Yes	Yes	No
LAL16B-2-7-18-06	Soil	High	28.94	2.91E+05	\pm	1.95E+04	Yes	Yes	No
RS-901-BJO	Water	Low	1.3	Not present			Yes	N/A	N/A
LAL16-GW2-7-18-06	Water	None	2.0	Not present			Yes	N/A	N/A
CollinsWell#1-7-7-06	Water	None	4.0	Not present			Yes	N/A	N/A
66783-7-26-06	Water	None	0.8	Not present			Yes	N/A	N/A
LK-01-0-01-8-9-06	Water	None	5.2	Not present			Yes	N/A	N/A
Hester-498-8-10-06	Water	None	2.9	Not present			Yes	N/A	N/A
N/A, not applicable. The sample was not run with the nested qPCR assay and/or the biomarker melt peak was not identified because none was detected in the qPCR sample run.									
Inhibited indicates that the sample did not amplify with qPCR even after a sepharose cleanup was performed and the sample was diluted to a lower DNA concentration.									
* Concentration in water is given in ng/L, concentration in soil is in ng/g									

4.4 Determine the Effect of Dilution on Quantification of the LA35 Biomarker in Soil and Water Samples by the qPCR Assay

4.4.1 Objective

This test was conducted to determine a concentration for poultry litter in an environmental sample below which the biomarker is likely not to be detected in environmental samples (i.e., simulate runoff conditions in the lab and test for the presence of biomarker at increasingly dilute samples).

4.4.2 Methods

Litter samples FAC2 and FAC8 (see Table 5) were each carried through the entire procedure. Note that FAC2 was one of the original litter samples used to find LA35. Litter sample FAC8 is a turkey litter sample. Both FAC2 and FAC8 previously tested positive for the poultry litter specific biomarker by the qPCR assay. An outline of the protocol is presented in Figure 4.

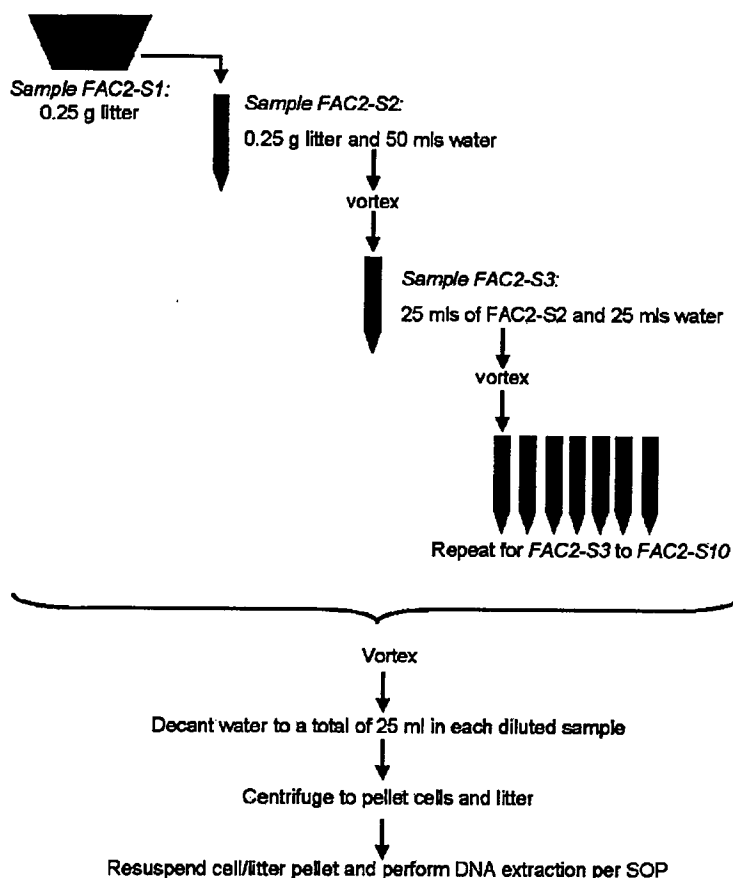


Figure 4. Outline of "dilution limit" protocol using litter sample FAC2 as the example.

4.4.3 Results

LA35 was amplifiable in samples down to 0.1 mg of litter per L. The concentration of LA35 in the litter samples themselves was greater than 2×10^8 copies/g of litter. There was a very strong correlation between litter concentration (grams of litter per L) and the concentration of LA35 (copies of LA35 per gram of litter), as indicated by the R^2 values of 0.97 and 0.99 shown in Figure 5. These R^2 values indicate that the DNA extraction is efficient and *Brevibacteria* quantification method is relatively precise. Table 9 shows the results of the analysis.

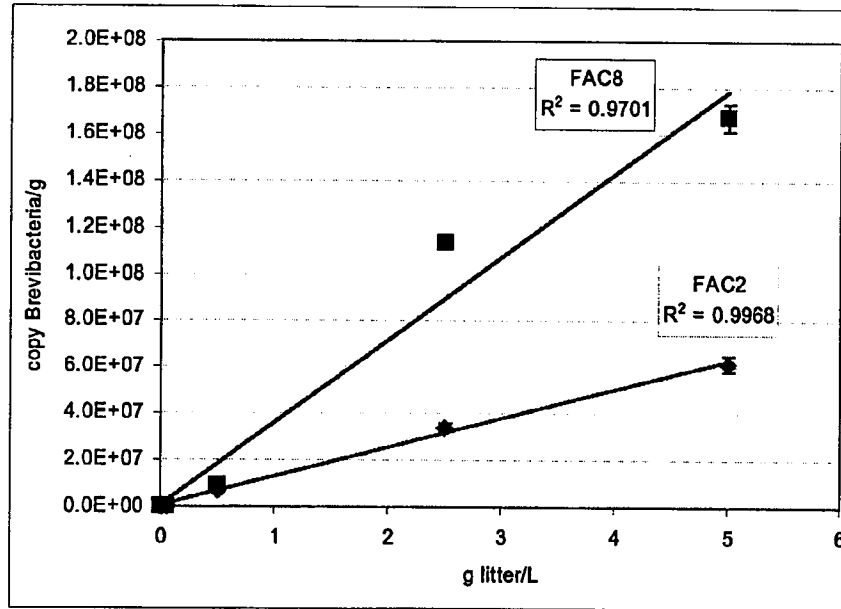


Figure 5. Copies of *Brevibacteria* poultry litter biomarker per gram of soil versus grams of soils per litter in the extractions.

Table 9. Results of the PCR, qPCR, and nested qPCR for *Brevibacteria* in litter.

Sample	"Litter" Concentration [^]	PCR Bacteria [#]	Nested qPCR [*]	<i>Brevibacteria</i> 16S rRNA (copy/L water or g litter)
FAC2-S1	NA	+	NA	2.32E+08 ± 8.88E+06
FAC2-S2	5,000 mg/L	+	NA	6.14E+07 ± 3.37E+06
FAC2-S3	2,500 mg/L	+	NA	3.42E+07 ± 1.61E+06
FAC2-S4	500 mg/L	+	NA	6.86E+06 ± 8.27E+05
FAC2-S5	50 mg/L	-	NA	3.54E+04 ± 4.84E+03
FAC2-S6	10 mg/L	-	NA	1.90E+04 ± 4.79E+03
FAC2-S7	5 mg/L	-	NA	1.01E+04 ± 9.35E+03
FAC2-S8	1 mg/L	-	NA	Present, not quantifiable

Table 9. (continued).

Sample	"Litter" Concentration [^]	PCR Bacteria [#]	Nested qPCR [*]	<i>Brevibacteria</i> 16S rRNA (copy/L water or g litter)
FAC2-S9	0.1 mg/L	-	+	Present, not quantifiable
FAC2-S10	0.01 mg/L	-	-	Not detected
FAC8-S1	NA	+	NA	2.56E+08 ± 2.49E+07
FAC8-S2	5,000 mg/L	+	NA	1.67E+08 ± 5.88E+06
FAC8-S3	2,500 mg/L	+	NA	1.14E+08 ± 1.56E+05
FAC8-S4	500 mg/L	+	NA	9.05E+06 ± 8.69E+05
FAC8-S5	50 mg/L	-	NA	7.48E+03 ± 2.53E+03
FAC8-S6	10 mg/L	-	NA	7.59E+04 ± 2.97E+04
FAC8-S7	5 mg/L	-	NA	3.20E+04 ± 2.23E+04
FAC8-S8	1 mg/L	-	+	Present, not quantifiable
FAC8-S9	0.1 mg/L	-	-	Not detected
FAC8-S10	0.01 mg/L	-	-	Not detected
[^] NA indicates that this was the original 0.25 g undiluted litter sample from which the dilutions were made. [#] "+" indicates that PCR products were observed after gel electrophoresis, indicating that amplifiable products were obtained from the sample. "-" indicates that PCR products were not observed by gel electrophoresis and UV visualization. Note that the detection limits of the gel electrophoresis method are much higher than the qPCR method detection limits. [*] "+" indicates that the <i>Brevibacteria</i> were detected with the nested qPCR protocol. "NA" indicates that the samples were not tested with the nested qPCR protocol.				

5 REFERENCES

- Bernhard, A.E. and K.G. Field, 2000, "Identification of nonpoint sources of fecal pollution in coastal waters by using host-specific 16S ribosomal DNA genetic markers from fecal anaerobes," *Applied and Environmental Microbiology*. 66(4): p. 1,587-1,594.
- Tsen, H, et al., 1998, "Development and use of 16S rRNA gene targeted PCR primers for the identification of *Escherichia coli* cells in water," *Journal of Applied Microbiology*. 85(3): p. 554-560.



Poultry-Specific Biomarker Quantitative PCR Analytical Summary

November 3, 2007

Overview:

The objective of this project was to quantify the number of poultry-specific *Brevibacteria* biomarker gene copies contained in water, soil, and/or litter samples using quantitative polymerase chain reaction (qPCR). Table 1 describes the sample matrix and the condition of the samples upon arrival to the analytical laboratory.

Table 1. Description of samples and volume or mass filtered for DNA extraction.

Sample ID	Matrix/ Date Sampled	Condition Received/Observations	Volume Filtered (mL) or Mass Extracted (g)
FAC03-7-6-06	Litter/7-6-06	Cold/sealed bag	0.35 g
FAC09-8-31-06	Litter/8-31-06	Cold/sealed bag	0.22 g
FAC02-6-21-06	Litter/6-21-06	Cold/sealed bag	0.24 g
FAC08-8-15-06	Litter/8-15-06	Cold/sealed bag	0.26 g
FAC04-7-12-06	Litter/7-12-06	Cold/sealed bag	0.26 g
FAC010-9-22-06	Litter/9-22-06	Cold/sealed bag	0.2 g
FAC07-8-3-06	Litter/8-3-06	Cold/sealed bag	0.39 g
FAC05-7-13-06	Litter/7-13-06	Cold/sealed bag	0.35 g
FAC06-7-20-06	Litter/7-20-06	Cold/sealed bag	0.27 g
LAL6A2-6-14-06	Soil/6-14-06	Cold/sealed bag	0.24 g
LAL6D2-6-15-06	Soil/6-15-06	Cold/sealed bag	0.24 g
LAL6B2-6-14-06	Soil/6-14-06	Cold/sealed bag	0.23 g
LAL13C2Q-7-7-06	Soil/7-7-06	Cold/sealed bag	0.24 g
LAL14C2-7-10-06	Soil/7-10-06	Cold/sealed bag	0.58 g
LAL14C2Q-7-10-06	Soil/7-10-06	Cold/sealed bag	0.51 g
LAL13C2-7-7-06	Soil/7-7-06	Cold/sealed bag	0.35 g
LAL16C2-7-18-06	Soil/7-18-06	Cold/sealed bag	0.35 g
LAL5C2-6-12-06	Soil/6-12-06	Cold/sealed bag	0.25 g
LAL8B2-6-21-06	Soil/6-21-06	Cold/sealed bag	0.26 g
LAL14B2-7-10-06	Soil/7-10-06	Cold/sealed bag	0.58 g
LAL9D2-6-22-06	Soil/6-22-06	Cold/sealed bag	0.25 g
LAL7B2-6-20-06	Soil/6-20-06	Cold/sealed bag	0.25 g
LAL13A2-7-6-06	Soil/7-6-06	Cold/sealed bag	0.27 g
LAL17C2-7-18-06	Soil/7-18-06	Cold/sealed bag	0.71 g
LAL17C2Q-7-18-06	Soil/7-18-06	Cold/sealed bag	0.62 g
LAL14D2-7-10-06	Soil/7-10-06	Cold/sealed bag	0.63 g
LAL8A2-6-19-06	Soil/6-19-06	Cold/sealed bag	0.7 g



Sample ID	Matrix/ Date Sampled	Condition Received/Observations	Volume Filtered (mL) or Mass Extracted (g)
LAL8D2-6-20-06	Soil/6-20-06	Cold/sealed bag	0.26 g
LAL17A2-7-10-06	Soil/7-10-06	Cold/sealed bag	0.68 g
LAL7D2-6-29-06	Soil/6-29-06	Cold/sealed bag	0.24 g
LAL9A2-6-22-06	Soil/6-22-06	Cold/sealed bag	0.26 g
LAL5A2-6-13-06	Soil/6-13-06	Cold/sealed bag	0.35 g
LAL7C2-6-19-06	Soil/6-19-06	Cold/sealed bag	0.25 g
LAL9B2-6-22-06	Soil/6-22-06	Cold/sealed bag	0.25 g
LAL13D2-7-6-06	Soil/7-6-06	Cold/sealed bag	0.24 g
LAL7A2-6-20-06	Soil/6-20-06	Cold/sealed bag	0.24 g
LAL16D2-7-18-06	Soil/7-18-06	Cold/sealed bag	0.52 g
LAL5B2-6-12-06	Soil/6-12-06	Cold/sealed bag	0.25 g
LAL12A2Q-7-6-06	Soil/7-6-06	Cold/sealed bag	0.27 g
LAL12D2-7-7-06	Soil/7-7-06	Cold/sealed bag	0.24 g
LAL16-SP2-7-18-06	Water/7-18-06	Cold/bottle intact	100 mL
EOF-SPREAD-010-5-9-06	Water/5-9-06	Cold/bottle intact	40 mL
EOF-SPREAD-17A-01-5-1-06	Water/5-1-06	Cold/bottle intact	30 mL
EOF-SPREAD-023-6-18-06	Water/6-18-06	Cold/bottle intact	25 mL
EOF-SPREAD-073B-6-18-06	Water/6-18-06	Cold/bottle intact	10 mL
EOF-SPREAD-064-5-4-06	Water/5-4-06	Cold/bottle intact	50 mL
EOF-SPREAD-53E-01-4-29-06	Water/4-29-06	Cold/bottle intact	30 mL
EOF-SPREAD-60-01-4-29-06	Water/4-29-06	Cold/bottle intact	50 mL
SPREAD-023-4-25-06	Water/4-25-06	Cold/bottle intact	40 mL
EOF-1-6-17-06	Water/6-17-06	Cold/bottle intact	100 mL
EOF-SPREAD-053G-5-4-06	Water/5-4-06	Cold/bottle intact	100 mL
EOF-SPREAD-048-5-9-06	Water/5-9-06	Cold/bottle intact	100 mL
SPREAD-029-4-25-06	Water/4-25-06	Cold/bottle intact	100 mL
SPREAD-036-4-25-06	Water/4-25-06	Cold/bottle intact	100 mL
EOF-SPREAD-071-5-9-06	Water/5-9-06	Cold/bottle intact	150 mL
EOF-SPREAD-065-5-4-06	Water/5-4-06	Cold/bottle intact	100 mL
EOF-Q2-6-17-06	Water/6-17-06	Cold/bottle intact	50 mL
EOF-26-6-8-05	Water/6-8-05	Cold/bottle intact	500 ml
EOF-17-6-8-05	Water/6-8-05	Cold/bottle intact	500 mL
EOF-222-4-13-07	Water/4-13-07	Cold/bottle intact	20 mL
GPGW-18A-6-26-07	Water/6-26-07	Cold/bottle intact	20 mL
GPGW-20-6-11-30-06	Water/11-30-06	Cold/bottle intact	50 mL
GPGW-48-7-12-1-06	Water/12-1-06	Cold/bottle intact	50 mL
GPGW-10-4-11-30-06	Water/11-30-06	Cold/bottle intact	50 mL



Sample ID	Matrix/ Date Sampled	Condition Received/Observations	Volume Filtered (mL) or Mass Extracted (g)
GPGW-40-6-27-07	Water/6-27-07	Cold/bottle intact	15 mL
HFS22-EVENTA-5-10-06	Water/5-10-06	Cold/bottle intact	150 mL
HFS04-BF2-01-8-1-06	Water/8-1-06	Cold/bottle intact	500 mL



Methods:

DNA Extraction. For soil and/or litter samples, DNA was extracted from 0.25 g of soil or litter using the FastDNA®SPIN® Kit for soil protocol. For surface water shipped to the laboratory, between 100 and 1,000 mL of groundwater was filtered through a Supor-200, 0.2 µm filter. The filters were frozen at -80°C and then shattered. Next, each sample tube was amended with 2 mL of DNA-free water, vortexed vigorously for 15 minutes, and the liquid volume was partitioned into DNA extraction tubes. DNA extractions were performed using the FastDNA®SPIN® Kit for soil according to the manufacturer's instructions. All DNA extractions were cleaned using an ethanol precipitation method. Community DNA was eluted in nuclease-free water (50 µL) and stored at -20°C.

Amplification of Bacteria. The PCR was used to amplify nearly full-length 16S rDNA genes from *Bacteria*. Each 25-µL PCR reaction included 1 X PCR buffer, 1.5 mM MgCl₂, 0.5 µM each 8F forward and 907R reverse primer, 1 u/50uL Taq DNA polymerase, 0.2 mM dNTP, 1 µL template DNA, and 20 µL molecular-grade water. Amplification was performed on a MJ Research Peltier Gradient thermocycler using the following regime: 94°C (5 min) followed by 30 cycles of 94°C (1 min), 53.5°C (1 min), and 72°C (1 min, 50 sec). The reaction was finished with an additional 7 minutes at 72°C. PCR products were examined by ultraviolet (UV) light in a 1% agarose gel stained with ethidium bromide to confirm specificity of the amplification reactions.

Sepharose Cleanup. Any sample not amplifying in the PCR was processed through a Sepharose CL-4B (Sigma-Aldrich) size exclusion gel chromatography cleanup. Briefly, the micro-bio spin columns (Bio-Rad) were packed with sterile Sepharose CL-4B and washed with Tris-HCl buffer (pH 8). The sample was added to the packed gel column and eluted by spinning in a micro-centrifuge.

Detection of a Poultry Specific *Brevibacterium* Biomarker. The qPCR methods for assessing the 16S rRNA gene are very sensitive in detecting specific DNA fragments. The detection limit for the methods used is approximately 6 gene copies per µL of the DNA extraction. Biomarker DNA was cloned into a plasmid and was used as the source of the quantitative standards used in the analysis. Plasmid DNA containing the target 16S rRNA gene from the poultry-specific *Brevibacterium* biomarker was purified and quantified fluorometrically. Based on the known size of the plasmid and insert, DNA concentrations were converted to insert copy numbers. A dilution series spanning seven orders of magnitude was generated using known concentrations of each plasmid. Amplification and detection of the DNA was performed using the MJ Chromo-4 System. The acceptance criterion for the standard curve is a linear R² value of greater than 0.995.

To determine qPCR results, sample DNA diluted to a final concentration of 15 ng/5 µL DNA was combined with following reagents to reach a final concentration of 1X SYBR Green Master Mix and 0.5 µM 157F and 727R primer and water to reach 20 µL and 5 µL, respectively, of diluted sample DNA. Amplification was performed on the MJ Research PTC-2004 thermocycler using the following regime: 50°C (2 min), 95°C (15 min), 40 cycles of 95°C (30 sec), 60°C (1 min), plate read and 50°C (5 min). The melting curve was determined using the following protocol: heat from 60°C to 90°C, by 0.3°C increments, holding for 5 seconds before reading the fluorescence of the samples. Nested qPCR results were determined by purifying the PCR products using the QIAquick PCR Purification Kit, as per the manufacturer's protocol, and then running the purified samples through qPCR, as described above.

QA/QC Requirements. To determine if and where potential contamination or interference occurred during sample processing, positives and reagent blanks or negatives and matrix spikes of the PCR and qPCR samples were prepared. A positive control consisting of pure DNA (known to amplify by specific DNA primers) was used for the PCR and qPCR procedure. A matrix spike consisting of pure DNA (known to amplify by specific DNA primers) was used for the PCR and qPCR procedure. Negative controls consisted of water-only blanks for the PCR and qPCR procedure. The qPCR reactions were run in triplicate for each sample to determine the reproducibility of the method.



Results:

The samples arrived at the lab in good condition at 4°C with ice still in the cooler. The samples were filtered in the lab, and the filters were immediately placed in a -80°C freezer and stored until the DNA extraction was performed. Following DNA extraction, the samples were first subjected to polymerase chain reaction (PCR) using universal bacterial probes in order to verify amplifiable DNA was present in the sample. In addition, for the 16S rRNA gene, a "nested" qPCR approach was applied in which the universal bacterial PCR-amplified DNA is used as the template in a qPCR reaction. Although the results from the nested qPCR cannot be quantified per se, they can be used to lower the detect limit for the qPCR in order to determine if the poultry-specific *Brevibacterium* biomarker gene is present at concentrations lower than the method detect limit (MDL) using the groundwater DNA extractions. The results of these studies are presented in Table 2. The DNA extraction negative control and all PCR negative controls did not amplify any product. In addition, all calibration control checks were within acceptable values.

Sample LAL16C2-7-18-06 results were reported previously (September 17, 2007). This sample was inhibited in the prior analysis as reported in September. This sample was sepharose cleaned an additional time and reanalyzed per the standard operating procedure (SOP). All inhibition was removed in the second cleanup and the *Brevibacterium* biomarker was detected in the qPCR analysis. The *Brevibacterium* biomarker was identified in 81% of the samples and was quantifiable in 39% of the samples analyzed.

Table 2. Results of Poultry Specific Biomarker analyses for samples.

Sample ID	Matrix	DNA (ng/L or ng/g)	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter) ^a	qPCR Matrix Spike Amplified? ^b	Nested qPCR Amplified? ^c	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
FAC03-7-6-06	Litter	21.3	1.03E+09 \pm 8.00E+07	Yes	N/A	Yes	No
FAC09-8-31-06	Litter	170.1	7.57E+08 \pm 1.55E+08	Yes	N/A	Yes	No
FAC02-6-21-06	Litter	51.9	4.13E+08 \pm 1.78E+07	Yes	N/A	Yes	No
FAC08-8-15-06	Litter	154.0	1.47E+09 \pm 2.25E+08	Yes	N/A	Yes	No
FAC04-7-12-06	Litter	6.8	1.67E+08 \pm 2.98E+07	Yes	N/A	Yes	No
FAC010-9-22-06	Litter	120.1	2.04E+09 \pm 4.14E+08	Yes	N/A	Yes	No
FAC07-8-3-06	Litter	98.1	2.49E+09 \pm 9.54E+07	Yes	N/A	Yes	No
FAC05-7-13-06	Litter	76.6	1.47E+09 \pm 1.93E+08	Yes	N/A	Yes	No
FAC06-7-20-06	Litter	57.1	4.46E+08 \pm 7.34E+07	Yes	N/A	Yes	No
LAL6A2-6-14-06	Soil	10.5	1.55E+04 \pm 2.57E+03	Yes	N/A	Yes	No



Sample ID	Matrix	DNA (ng/L or ng/g)	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter) ^a	qPCR Matrix Spike Amplified? ^b	Nested qPCR Amplified? ^c	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
LAL6D2-6-15-06	Soil	2.1	4.98E+03 \pm 1.88E+02	Yes	N/A	Yes	No
LAL6B2-6-14-06	Soil	12.4	Present	Yes	Yes	Yes	No
LAL13C2Q-7-7-06	Soil	27.3	Present	Yes	Yes	Yes	No
LAL14C2-7-10-06	Soil	20.3	Present	Yes	Yes	Yes	No
LAL14C2Q-7-10-06	Soil	77.3	Present	Yes	Yes	Yes	No
LAL13C2-7-7-06	Soil	3.4	Present	Yes	Yes	Yes	No
LAL16C2-7-18-06	Soil	9.5	1.42E+04 \pm 1.97E+03	Yes	N/A	Yes	No
LAL5C2-6-12-06	Soil	40.3	Present	Yes	Yes	Yes	No
LAL8B2-6-21-06	Soil	27.7	0.00E+00	Yes	No	N/A	N/A
LAL14B2-7-10-06	Soil	91.7	Present	Yes	Yes	Yes	No
LAL9D2-6-22-06	Soil	13.0	Present	Yes	Yes	Yes	No
LAL7B2-6-20-06	Soil	47.4	Present	Yes	Yes	Yes	No
LAL13A2-7-6-06	Soil	38.6	Present	Yes	Yes	Yes	No
LAL17C2-7-18-06	Soil	17.3	Present	Yes	Yes	Yes	No
LAL17C2Q-7-18-06	Soil	29.9	Present	Yes	Yes	Yes	No
LAL14D2-7-10-06	Soil	28.1	Present	Yes	Yes	Yes	No
LAL8A2-6-19-06	Soil	22.3	7.00E+03 \pm 4.43E+02	Yes	N/A	Yes	No
LAL8D2-6-20-06	Soil	38.4	0.00E+00	Yes	No	N/A	N/A
LAL17A2-7-10-06	Soil	21.4	Present	Yes	Yes	Yes	No
LAL7D2-6-29-06	Soil	20.6	Present	Yes	Yes	Yes	No
LAL9A2-6-22-06	Soil	29.2	Present	Yes	Yes	Yes	No
LAL5A2-6-13-06	Soil	35.9	Present	Yes	Yes	Yes	No
LAL7C2-6-19-06	Soil	41.8	Present	Yes	Yes	Yes	No
LAL9B2-6-22-06	Soil	61.4	Present	Yes	Yes	Yes	No
LAL13D2-7-6-06	Soil	34.6	Present	Yes	Yes	Yes	No



Sample ID	Matrix	DNA (ng/L or ng/g)	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter) ^a	qPCR Matrix Spike Amplified ^b	Nested qPCR Amplified ^c	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
LAL7A2-6-20-06	Soil	31.7	Present	Yes	Yes	Yes	No
LAL16D2-7-18-06	Soil	36.4	Present	Yes	Yes	Yes	No
LAL5B2-6-12-06	Soil	27.3	Present	Yes	Yes	Yes	No
LAL12A2Q-7-6-06	Soil	26.8	Present	Yes	Yes	Yes	No
LAL12D2-7-7-06	Soil	51.4	Present	Yes	Yes	Yes	No
LAL16-SP2-7-18-06	Water	-1.0	0.00E+00	Yes	No	N/A	N/A
EOF-SPREAD-010-5-9-06	Water	1.7	1.05E+07	Yes	N/A	Yes	No
EOF-SPREAD-17A-01-5-1-06	Water	72.5	2.48E+06	Yes	N/A	Yes	Yes
EOF-SPREAD-023-6-18-06	Water	4.3	1.11E+05	Yes	N/A	Yes	No
EOF-SPREAD-073B-6-18-06	Water	133.5	5.56E+07	Yes	N/A	Yes	No
EOF-SPREAD-064-5-4-06	Water	79.2	1.89E+06	Yes	N/A	Yes	No
EOF-SPREAD-53E-01-4-29-06	Water	57.7	5.45E+07	Yes	N/A	Yes	No
EOF-SPREAD-60-01-4-29-06	Water	431.4	3.90E+07	Yes	N/A	Yes	No
SPREAD-023-4-25-06	Water	194.2	1.25E+06	Yes	N/A	Yes	No
EOF-1-6-17-06	Water	2.5	1.15E+05	Yes	N/A	Yes	No
EOF-SPREAD-053G-5-4-06	Water	14.5	0.00E+00	Yes	No	N/A	N/A
EOF-SPREAD-048-5-9-06	Water	25.2	0.00E+00	Yes	No	N/A	N/A
SPREAD-029-4-25-06	Water	56.1	Present	Yes	Yes	Yes	No
SPREAD-036-4-25-06	Water	64.9	1.48E+05	Yes	N/A	Yes	No
EOF-SPREAD-071-5-9-06	Water	5.2	3.63E+04	Yes	N/A	Yes	No
EOF-SPREAD-065-5-4-06	Water	13.3	3.45E+04	Yes	N/A	Yes	No
EOF-Q2-6-17-06	Water	2.0	0.00E+00	Yes	No	N/A	N/A
EOF-26-6-8-05	Water	14.3	0.00E+00	Yes	No	N/A	N/A
EOF-17-6-8-05	Water	22.7	Present	Yes	Yes	Yes	No
EOF-222-4-13-07	Water	81.3	1.32E+05	Yes	N/A	Yes	No
GPW-18A-6-26-07	Water	1.1	0.00E+00	Yes	No	N/A	N/A



Sample ID	Matrix	DNA (ng/L or ng/g)	qPCR Poultry Specific Biomarker (copies/ μ L water or g soil or g litter) ^a	qPCR Matrix Spike Amplified? ^b	Nested qPCR Amplified? ^c	Biomarker Melt Peak Identified?	Other Melt Peaks Observed?
GPGW-20-6-11-30-06	Water	13.4	0.00E+00	Yes	No	N/A	N/A
GPGW-48-7-12-1-06	Water	2.6	0.00E+00	Yes	No	N/A	N/A
GPGW-10-4-11-30-06	Water	2.9	Present	Yes	Yes	Yes	No
GPGW-40-6-27-07	Water	-1.5	0.00E+00	Yes	No	N/A	N/A
HFS22-EVENTA-5-10-06	Water	4.2	0.00E+00	Yes	No	N/A	N/A
HFS04-BF2-01-8-1-06	Water	0.0	0.00E+00	Yes	No	N/A	N/A

^a "Present" indicates that the biomarker was amplified, but was not quantifiable.
^b If "no" indicates that sample did not amplify with qPCR even after a sepharose cleanup was performed and the sample was diluted to a lower DNA concentration indicative of inhibition.
^c N/A, not applicable. The sample was not run with the nested qPCR assay and/or the biomarker melt peak was not identified because none was detected in the qPCR sample run.